



FlexPoint® Piercing Stylet (PS)

The Piercing Stylet is an accessory, used in conjunction with a matched size STSN device, to increase the tip sharpness for crossing tissue that is unusually scarred and difficult to penetrate.


Advantages

- Leverages current technique
- Inserts within the lumen of the STSN before or during the procedure
- Locks to proximal end of STSN
- Safety Clip prevents inadvertent activation until required
- Spring-loaded system retracts immediately upon release of activation button
- Does not leave permanent RF/cautery damage to the septum
- Can be quickly removed from the STSN for pressure measurement and dye injection


Available Configurations

Description	Working Length	Model Number	Part Number
Piercing Stylet for 40-cm STSN	49.5-cm	PS-40	FA07-0017
Piercing Stylet for 71-cm STSN	80.5-cm	PS-71	FA07-0018
Piercing Stylet for 89-cm STSN	98.5-cm	PS-89	FA07-0019
Piercing Stylet for 98-cm STSN	107.5-cm	PS-98	FA07-0020


PS-71 Projecting Beyond STSN Tip

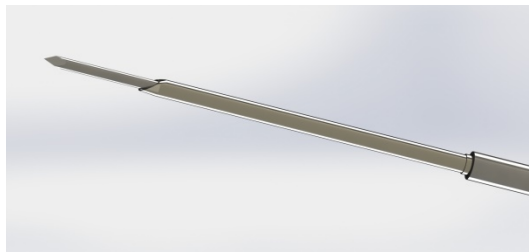

PS-71 Hub with Safety Clip

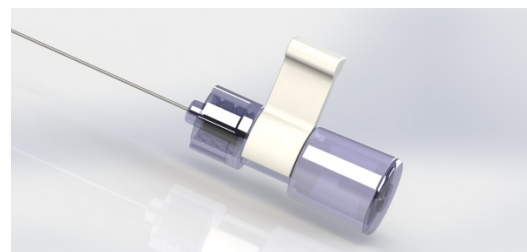

PS-71 on STSN with Safety Clip Removed

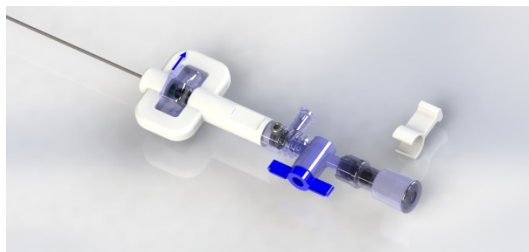

PS-71 on STSN Hub Activated

